

BEAUTECH



BEAUTY DEVICE

Read this manual carefully before use & keep it for future reference.

Products

Thank you for purchasing this RF Beauty Device.

This product is a revolutionary new radio frequency beauty instrument. The use of RF technology in skincare for multi-dimensional anti-aging focuses on rejuvenating deeper skin layers, offering a more holistic and effective solution for reducing the visible signs of aging.

The stamping-type radiofrequency technology represents a breakthrough in RF treatments by improving precision, uniformity, and depth of penetration. It effectively targets the dermis layer to stimulate collagen production, tighten the skin and reduce the visible signs of aging, offering a powerful, non-invasive solution for anti-aging skin rejuvenation.

The energy generated by traditional radiofrequency instruments often only acts on the surface layer of the skin, making it difficult to act on the dermis. This product with radio frequency technology with stamping type of approach can make the radio frequency energy more focused, penetrating into the dermis to significantly increase the amount of collagen regeneration, so as to effectively achieve lifting and tightening, lifting the jawline, reduce the effectiveness of wrinkles. In addition, to the stamping mode, the traditional sliding mode is also retained to meet the usage habits of different needs.

Multi-mode technology, RF technology and EMS technology, and micro vibration massage at the same time, helps the skin better absorb the nutrients of skincare products, and lift and tighten the skin from the inside out to get better skincare results, so that your skin always show a healthy, moisturised and elasticity.

Caveat

Please read this user manual carefully and keep it in a safe place for future reference.

Important

- Do not allow children to touch this product.
- Do not throw the instrument into fire.
- Keep the product in a safe place to prevent it from being dropped and damaged.

Warnings

- When charging, only use a charger with 5V output or a mobile power supply to charge this instrument.
- Do not continue to use this instrument if it appears damaged to avoid injury.
- Completely charge the instrument at least once every 2 months to protect the battery life.
- This instrument is not intended to be used by children under 12 years of age or individuals with physical, sensory, or cognitive impairments.
- Do not insert metal-containing materials into the jack of the instrument to avoid short-circuiting.
- It is recommended that women who are pregnant or breastfeeding should not use this instrument.

Caution


- Do not use while sleeping or sleepy.
- Do not use this instrument for any other purpose than that described in this manual.
- It is recommended to charge, use and store this instrument at temperatures between 5°C/40°F and 35°C/95°F.
- Do not use the device if you are taking steroid medication.
- Do not use this device if your skin is still recovering from a cosmetic procedure (e.g., chemical or mechanical exfoliation, laser resurfacing, fractional laser treatments, electrolysis treatments, or any other filler-contouring treatments performed within the last 2 weeks)

Warnings

- If you have any medical concerns, consult a dermatologist first before use.
- Do not use this device on severe nodular or cystic acne.
- Do not use the device on moles, warts, ingrown hairs or irritated skin.
- Do not use the device if you are sensitive to light.
- Do not use this device if you have had sunbeam or skin burns in the past two weeks.
- Do not use the device over the eyeballs.
- Do not use the same area more than 3 times in a day. Excessive use may cause burns or damage to the skin. Allow at least 4 hours between treatments.
- If you experience any skin discomfort (tightness, discomfort, irritation, redness, etc.) after using the device, please do not use the

- device on the same area in a day. (comfort, irritation, redness, burning, etc.), notice unexpected side effects or prolonged redness, stop using it immediately and consult your doctor. However, temporary slight redness is normal. However, temporary slight redness is normal.
- Do not look directly at the light source of the probe for a long time, as this may cause eye damage.
- Do not allow children to play with this device.
- Do not share your instrument with others as this may result in cross contamination of bacteria.
- If the instrument appears visibly damaged, refrain from use and avoid attempting to disassemble or repair it yourself.

product mix

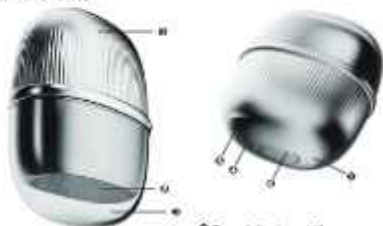

- ①Hosts
- ②Dot matrix probe
- ③Dust cover

- ④Gear Adjustment Keys
- ⑤Power on key/mode switching key
- ⑥Charging port
- ⑦Gear indicator/charging indicator

Functional guidelines

1. The instrument has three modes: Stamping mode, EMS mode and Gold Slide Lift mode.
 2. Long press the power on button to turn on the instrument, you can select the desired function mode by mode switching.
 3. After the instrument is turned on, you can select the intensity that suits you through the gear switching key.
 4. Long press the on/off button for 3 seconds to turn on or off the instrument.
- *Automatic protection shut-off function When you do not have any key operation, the instrument will turn off automatically after 5 minutes of continuous use.
5. To turn on the vibration double-click the on/off button; if you double-click the on/off button again, the vibration will be cancelled; when switching modes, the vibration will also be turned off automatically.

Stamping mode

1. Principle and efficacy The RF radio frequency will work simultaneously with nano red light.

Radiofrequency (RF) technology works by delivering RF waves into the subcutaneous tissue, where the natural resistance of the tissue to heat causes a contraction of collagen fibres. This process tightens the skin while simultaneously promoting collagen regeneration. The immediate contraction of collagen fibers helps to lift and firm the skin, while the stimulation of the dermal layer encourages the production of new collagen, improving the thickness and density of the skin's dermal layer.

This increased dermal thickness and density reduce wrinkles and improve skin laxity.



Additionally, using LED nano red light before or after the treatment enhances the microcirculation of the lymphatic system and stimulates mitochondrial activity within cells. The specific wavelength of the red light neutralises free radicals and encourages fibroblasts to produce more collagen. This helps to accelerate blood circulation and provides photo-rejuvenation effects.

The continuous vibration wave produced by the device also massages the skin, promoting nutrient absorption, tightening the skin, and effectively preventing wrinkles and sagging.

2.How to use

Step 1: Long press the On button to turn on the device, the stamping mode will start by default, and the default gear is 1st. Adjust the gear according to the degree of stimulation your skin will receive. The higher the gear, the better the effect.

Step 2: After applying nutrient essence (or when applying mask), use the massage probe to touch the skin to perform a stamping and pressing operation.

Step 3: For the face in accordance with the sequence of stamping, press operation, after 8 seconds, there will be a "ticking" sound, indicating that the stamping is complete, change to the next stamping region, until the full face is complete.

Description: The radiofrequency wave generated by this device features a built-in safety mechanism. When the tone (a "beep") sounds, the radiofrequency energy automatically decreases, and the device enters a semi-dormant protective state. If the device is moved away from the skin after the "beep", the energy will remain in this safe state. When the device is re-applied to the skin for continued treatment, the radiofrequency energy is activated again. This design ensures that even if the device is not removed from the skin immediately after the "beep" reminder, it will not cause harm to the skin.

EMS model

1.Principle and efficacy EMS microcurrent will work simultaneously with nano blue cosmetic light.

EMS microcurrent can revitalise facial muscles and epiderma nerves.

Improve movement, activate cells, stimulate cell regeneration, promote metabolism, slow down skin puffiness, lifting and firming the skin.

LED nano blue light adopts 465nm wavelength, which can regulate the oil-water balance of the skin, and it has the effect of cleaning and sterilising, purifying skin and eliminating acne.

2.How to use

Step 1: Long press the on button to turn on the machine, it will default to stamping mode (red light), tap the on/off button to switch to blue light which is the EMS mode, the default gear is 1. The higher the gear, the more stimulation to the skin.

Step 2: After applying the essence (or when applying mask), touch the skin with the massage probe in a slow sliding operation.

Operation Recommendations

Gear recommendation:

①Thin skin in the forehead area, it is recommended to use grades 1-2

③Moderate thickness in the cheek area, 2-3 speeds recommended

②The skin around the eye area is the thinnest, it is recommended to use grades 1-2

④Moderate skin thickness in the jaw area, 2-3 speeds recommended



*It is recommended that [Stamp Mode] be used for 5 minutes to promote collagen regeneration, EMS Model Use with the remaining essence or gel of the mask for 5 minutes to lift and firm.



Gold Slide Lift Mode

1. Principle and Efficiency

Radio Frequency will work simultaneously with Nano Amber Cosmetic Light.

RF radiofrequency waves generate heat energy to increase the thickness and density of the dermis layer of the skin, repairing reticular fibers, fading wrinkles, and improving elasticity.



LED Nano Amber Light uses 588-598nm wavelength, which can improve the skin's microcirculation and increase skin elasticity, soothe wrinkles, stretch marks, and slow down skin aging.

WARNING: The RF waves generated by the RF radio frequency in this mode will continuously heat the dermal layer of the skin, please keep gliding the massage probe when it is attached to the skin to operate, not too much to stay, stay for a maximum of 6 seconds, do not stay in the same area of the skin for too long, so as not to burn the skin.

2. How to use

Step 1: Long press the power on button to turn on the machine, the default first start stamping mode (red light), tap the power on button, that is, switch to [sliding mode] (amber light), the default is 1 gear, you can according to their own skin to receive stimulation, gear adjustment; the larger the gear the better the effect.

Step 2: After applying the nutritional essence (or when applying the mask), use the massage probe to contact the skin for slow sliding operation.

[Stamping mode]:

It is recommended to use 1-4 times a week depending on individual skin type.

- Normal skin - 10-15 minutes 3 times a week
- Dry skin - 10-15 minutes 2 times a week
- Oily skin - 15-20 minutes 3-4 times a week
- Sensitive skin - 5-10 minutes 1 time per week

[EMS model]:

Recommended for use every 2 days.

[Stamped Anti-Aging Mode] and [EMS Skin Contouring Mode] can be used in conjunction with each other for a complete beauty treatment. It is recommended to use [Stamped Anti-Aging Mode] for 5 minutes to promote collagen regeneration and [EMS Contouring Mode] for 5 minutes to lift and firm with the remaining essence or gel of the mask.

[Gold Slide Lift Mode]:

It is recommended to use it once every 2 days.

[Golden Slide Lifting Mode] This mode is a mixed mode of RF+EMS, which can be used to lift the jaw, and can be used selectively according to personalised needs.

The recommended order of using the facial areas for all three modes is: left cheek, right cheek, and forehead.

Treatment step

Step 1: Remove Makeup and thoroughly cleanse your face and neck to ensure your skin is in a clean condition.



Step 2: Apply a Gel or Mask

[Stamped Anti-Aging Model]:

Remove the mask, unfold the mask and apply it to the area you want to use on the face (the mask should sit flat to the face).

Or apply essence on the face, or you can directly squeeze the essence (gel) directly onto the instrument probe.



Step 3: Facial Treatment

Power on the unit and select the desired mode to either use on its own or switch between the 3 modes one at a time.

When there is no key operation during use the instrument turns off automatically, indicating that the mode has been used for 5 minutes. Please reboot the instrument and continue to use it.



Step 4: Wash and hydrate the skin

At the end of the treatment, wash your face with water and apply the right amount of skin care products to hydrate and moisturise according to your personal needs.

Step 5: Machine maintenance

- After the treatment, please use a slightly dampened soft cloth to wipe off the residue on the beauty instrument.
- If there is a lot of residue, you can rinse it directly with water and after rinsing, dry it with a dry cloth.



Warning: Although this product can be rinsed, over-rinsing may reduce the service life of the equipment, please try to wipe clean.

How to Charge

1. Connect the USB charging cable to a charging head with a 5V output (or to a computer or portable power source) and plug the Type-C connector into the charging cradle port.
2. Once the device is fully charged, the indicator light will become solid and stop pulsing or undulating. This indicates that the charging process is complete, and the device is no longer drawing power, as the battery is full.

Attention

Before charging, make sure that the charging port is dry and free of any residing skin care product. Do not turn on while the device is charging.

Cleaning and maintenance of instruments

- Do not use abrasive cleaners or caustic liquid (e.g. Acetone) to clean the instrument.
- Thoroughly clean the massage probe with water after each use and allow to dry naturally to ensure maximum performance.

Table of Common Failure Problems

concern	Presumed cause	Solutions
Unable to start the device	Battery depletion	Use after charging is complete
	Battery life depletion	Contact seller for repair or replacement
Automatic shutdown of equipment	No operation for more than 5 minutes	Restart
	Battery depletion	Use after charging is complete
	Not enough essence	Replenish the amount of essence on your face
Unable to charge	Use of mismatched charging devices	Use a charging plug with 5V output or rechargeable battery)
Feeling that your face is too hot	Gear is too high	Reduced gear strength
	The beauty device stays in the area and does not move	Slightly faster product movement

Feeling that your face is not hot enough	Too low a gear	Boosted gear strength
	Groomer slides too fast	Slightly reduces product movement speed
	Not enough essence	Regulate the amount of essence on your face
	Insufficient energy storage when the unit is first turned on	Let the product run a little longer
Tingling current	EMS gears are too strong	Turn down a gear
	Device is not skin-tight	Holding the device close to the skin
Tingling sensation	Only the edge of the probe on the skin	Probe flat against the skin
Slight redness of the skin after use	Normal phenomenon of RF wave action	Normalization can be achieved in a short period of time, if there is non-recovery or other adverse reactions. Please consult a doctor promptly
Water always flows out during use	Rinse products without water control	Pat the product gently on a dry cloth to control the moisture or let it dry in a ventilated area before use
Beeping sounding use	Different modes	Only the stamping mode has a "beep" sound, the other modes do not have a sound, please read the manual in detail

Statement denying or limiting responsibility

Product users do so at their own risk. No responsibility or liability is assumed for any physical or non-physical injury resulting directly or indirectly from the use of the product. In addition, the company reserves the right to modify these instructions and change the contents at any time without any obligation to notify any person of such modifications or changes.

Disposal methods

Waste disposal methods for used electronic equipment (applicable in the EU or countries or regions with waste separation systems).



The pattern of forks painted on the disposal bin indicates that the instrument should not be disposed of as household waste and should be disposed of in a suitable electrical and electronic instrument recycling collection station (or recycling bin). Proper disposal of this instrument prevents waste instruments from negatively impacting the environment and human health. Material recycling also helps to conserve natural resources. It is everyone's responsibility to protect the Earth.

Product Specification

Product Name: BEAUTY DEVICE

Product powers: 10W

Frequency: 1000KHz(1MHz)

Battery capacity: 1100mA

Charging time: 3H

Input Voltage: 5V

USB charging interface: Type-C

Waterproof rating: IPX6

Theoretical usage time for a full charge:

Stamping mode: about 60 minutes

EMS mode: approx. 60 minutes

Gold Slide Lift mode: approx. 120 minutes

*Depending on the selected gear and different skin resistance value, it will vary

Standby time: 60 days

NOTE

- Substance of the following statement: "This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety."
- Do not use this heater in the immediate surroundings of a bath, a shower or a swimming pool.
- Do not use this heater if it has been dropped
- Do not use if there are visible signs of damage to the heater use this heater on a horizontal and stable surface.
- WEEE disposal information



Worked 4 Subdivisions
880P 101 110
Phone: 1800-111-1111
Address: 1000000, 1000000, 1000000
Tel: 10000000000
Email: 10000000000@1000000.com

MADE IN CHINA

