


BEAUTECH

## RF EMS Therapy Device User Manual

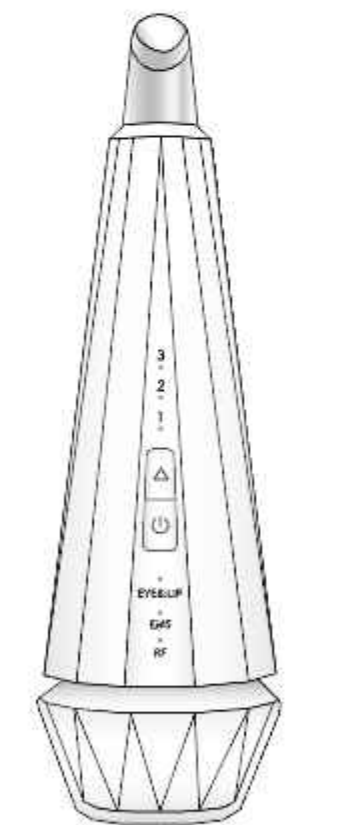

Thank you very much for purchasing this product  
please read the manual carefully before use


### Manual

Product Specifications	1
Product Structures	1
Packaging and Accessories	2
Product Features	2
Functions	2
Mode Introduction	4
How to Use	4
Caution	6
Warning	7


### Product Specifications

Model No.: KFB259  
Charging Type: USB Charging  
Power Voltage: 3.7V  
Product Power: 2.59W  
Battery Capacity: 1000mAh  
Charging Duration: 2.5 hours  
Working Times: total 2 hours  
Standby Time: 90 Days (If the unit hasn't been used for more than 90 days, fully charged before use)  
Material: PC + Zinc Alloy  
Vibration Frequency: 625Hz  
Power Supply: DC 5V 1A

### Product Structures

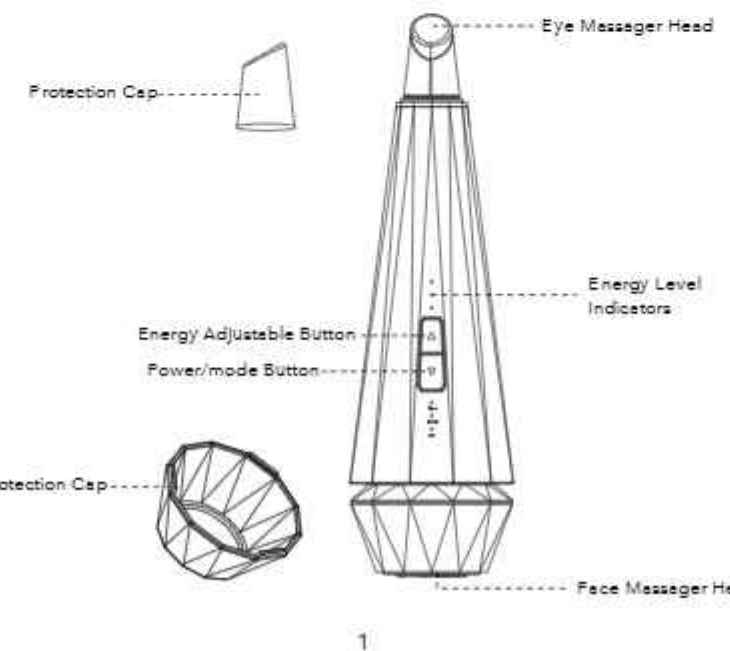

### Back

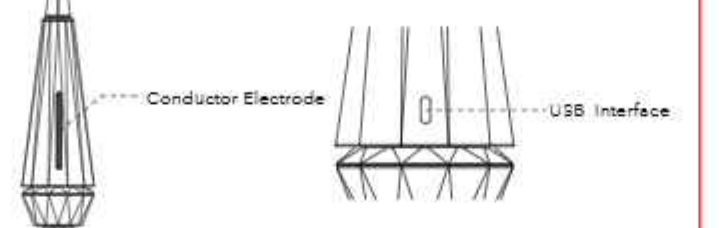

### Packaging and Accessories

1. Device
2. USB Cable
3. User Manual
4. Cotton Pads

### Product Features

1. Dual-solution: The dual-solution device pairs a face massage head with eye massage head.
2. Diamond-shaped design for a luxurious and premium look.
3. Heating therapy directly relaxes muscles, improves blood circulation, and stimulates it.
4. Eco-friendly built-in rechargeable lithium battery.
5. The device features a high-quality metal massage head, durable polycarbonate (PC) material, and a solid structure.

### Functions

#### 1. RF (Radio Frequency)

Radio Frequency operates by delivering energy through closely spaced electrodes, allowing efficient penetration into the skin. When the high frequency is flowing between the tissues of living body, energy is produced in the easily flowing portion according to the characteristics of high frequency and the resistance of biological tissue.

#### 2. EMS (Electrical Muscle Stimulation)

EMS is used to activate the skin cells, stimulating collagen rejuvenation, to improve skin elasticity, soften wrinkles and leave skin looking smooth and healthy. Unique EMS current can activate cells, promote the proliferation of collagen protein, and stimulate muscles to lift and tighten, making the skin feel firm and elastic. Stimulation of subcutaneous tissue activates the cells and stimulate collagen to contract and reorganise, which create new collagen, making facial muscles firmer, more dynamic, and elastic. Therefore, this device can improve fine lines and wrinkles on the skin to regain a youthful and smooth appearance, without any damage to skin.

#### 3. Vibration

Vibrations acting on the face can revitalise skin cells through gentle vibration and create a micro-massage effect, increasing cell volume. This improves local blood and lymph circulation, strengthens cellular permeability, improves tissue metabolism and regeneration, softening tissue, stimulates the nervous system and cell function, leaving the skin radiant and elastic.

#### 4. Heating Therapy

Also called thermotherapy, is the use of heat in therapy, for pain relief and health. Heat applied to the skin directly opens pores. By opening the pores of the skin, it helps absorb better treatments such as facial and body masks. It also reduces joint stiffness; relieves pain; relieves muscle spasm; reduces inflammation, assists therapy acute edema, and increases blood flow.

#### 5. Negative Ion

Negative ionisation treatments which then bind to positive ions and the free radicals in the skin cells, which prevent cellular damage. This boosts collagen production and blood circulation.

### Mode Introduction

The device will turn off automatically at any mode in 10 minutes. Press Power / Mode Button to adjust mode.

#### EYE & LIP Mode (Eye & Lip Massage Head):

Heating, EMS, Negative Ion and Vibration will work together, which is fit for lip and eye area.

#### EMS Mode (Face Massage Head):

EMS, Negative Ion and Vibration work together, making it suitable for the face, forehead, and neck.

#### RF Mode (Face Massage Head):

RF, Negative Ion and Vibration are working making it suitable for the face, forehead and neck.

Note: There is only one level for EYE & LIP Mode, and 3 adjustable levels for EMS Mode and RF Mode.

### How to Use

#### EYE & LIP Mode

1. After makeup removal, apply toner, facial cream, serum, or eye cream evenly over the skin.
2. Please take off the Protection Cap of the Eye & Lip Head, then hold for two seconds to activate the default Eye & Lip mode, while white LED will be on of RF Mode and the first Energy Level Indicator.
3. Put the Massage Head on the skin and move it slowly as per the direction shown in the right picture.
4. After use, press the Power/Mode button to turn off the device and clean it with cloth.

#### EMS Mode

1. After makeup removal, apply toner, facial cream, serum, or eye cream evenly over the skin.
2. Please take off the Protection Cap of the Eye & Lip Head, then press the Power Button and press again to enter EMS mode, while white LED will be on of EMS Mode and the first Energy Level Indicator. Press Energy Adjustable Button to Adjust the levels.
3. Place the Massage Head on the skin and move it slowly in the direction shown in the picture on the right.
4. After use, press Power/Mode button to turn off the device and clean it with cloth.

#### RF Mode

1. After makeup removal, apply toner, facial cream, serum, or eye cream evenly over the skin.
2. Please take off the Protection Cap of the Eye & Lip Head, then press the Power Button and press twice to enter RF mode, while white LED will be on of RF Mode and the first Energy Level Indicator. Press Energy Adjustable Button to Adjust the levels.
3. Place the Massage Head on the skin and move it slowly in the direction shown in the picture on the right.
4. After use, press Power/Mode button to turn off the device and clean it with cloth.

### How to Charge

- Ensure the device is fully charged before its first use. If the rechargeable battery is empty when you want to use the device, charge it for at least 15-25 minutes to be able to perform one treatment.
1. Make sure the appliance is switched off.
  2. Connect the USB cable to the device and adaptor, then insert the adaptor into the socket (The device cannot be used while charging. During charging, the mode LED light will flash).
  3. When the battery is fully charged, the mode LED will flash during charging.
  4. When the battery is low, the mode LED will flash four times and then turn off.

### Caution

1. Ensure the device is switched off before cleaning.
2. Please don't immerse it in water or clean with fiercely running water. Keep this device dry and away from moisture.
3. Do not use diluents, solvents, or other corrosive industrial products to clean the sonic head and hot massage head, as this could cause damage.
4. The treatment may lead to mild redness directly after the treatment, because the blood flow will be influenced. You may feel a warm, tingling sensation during the treatment. These effects usually disappear within several minutes and are generally perceived as pleasant.

### Warning

1. This device can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the device in a safe way and understand the hazards involved.
2. Children should not play with the device.
3. Cleaning and maintenance must not be performed by children without supervision.
4. Do not insert any material into the USB socket to avoid short circuit.
5. Do not use the device if the treatment head is broken, cracked, damaged or dirty. Always inspect the treatment head before use.
6. Do not apply pressure on your eyeball.
7. Pregnant women and people with any kind of heart condition should avoid this treatment. Discuss any allergy or skin sensitivities with your cosmetic doctor or aesthetician before starting a treatment.
8. This symbol means the device must not be cleaned under a running tap.
9. Keep the device and USB cable dry.
10. Substance of the following statements: This appliance is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
11. Do not use this device in the immediate surroundings of a bath, a shower, or a swimming pool.
12. Do not use this device if it has been dropped.
13. Do not use this device if there are visible signs of damage.
14. Use this device on a horizontal and stable surface.
15. Refer to WEEE disposal guidelines for proper disposal instructions.



Imported & Distributed by  
**Shop TJC Ltd.**  
Surrey House, Plane Tree Crescent  
Feltham, Middlesex, UK, TW13 7HF  
Tel: 08002606666  
Email: [contactus@tjc.co.uk](mailto:contactus@tjc.co.uk)