

BEAUTECH



BEAUTY DEVICE

Read this manual carefully before use & keep it for future reference.

Product Description

Thank you for purchasing this radiofrequency beauty device. This product is a revolutionary new radiofrequency (RF) beauty Device.

This product utilises advanced RF (Radio Frequency) technology to target both the epidermis and dermis for multi-dimensional anti-aging benefits. It offers a larger treatment area and enhanced levels of effectiveness, marking a significant breakthrough over traditional methods. By employing a stamped application technique, it delivers more focused, uniform intensity, ensuring deeper penetration of energy into the skin for improved anti-aging results.

The Stamped RF technology is a revolutionary approach where collagen renewal is achieved through heating the dermis to contract and tighten collagen cells beneath the skin. Traditional RF devices often only affect the surface layer, limiting their ability to reach the dermis. However, with the Stamped Mode, the RF energy is concentrated and penetrates deeper, significantly boosting collagen regeneration. This results in visible lifting and tightening effects, particularly along the jawline, and helps reduce the appearance of wrinkles.

Equipped with multi-mode technology, this device combines RF with EMS (Electrical Muscle Stimulation) and micro-vibration massage to enhance the absorption of skincare products and firm the skin. The combination lifts and tones the skin from the

inside out, leaving you with a healthy, radiant complexion. In addition to the Stamped Mode, the traditional sliding mode is also included, catering to the preferences of different users.

Warning

- When charging, only use a charger with an output of 5V or a mobile power supply to charge this instrument.
- Do not continue to use the instrument if it appears to be damaged as this may cause injury.
- Charge the apparatus fully at least once every 2 months to protect the battery life.
- This device is not suitable for use by children under 12 years of age, individuals with physical impairments, sensory or mental disabilities, or those lacking relevant experience or knowledge, unless they are supervised or instructed on its use by a responsible person. It is essential to supervise or provide instruction to ensure safe usage of the device.
- Do not insert materials containing metal into the jacks of the apparatus as this may cause a short circuit.
- Pregnant or breastfeeding women should avoid using this device.


Note

- Do not use while sleeping or sleepy.
- Do not use the instrument for any other purpose than that described in this manual.
- It is recommended to charge, use and store this instrument at temperatures between 5°C / 40°F and 35°C / 95°F.
- Do not use the device if you are taking steroids.
- Do not use the device if the skin is still recovering from a cosmetic procedure (e.g. chemical or mechanical exfoliation, laser resurfacing, fractional laser treatment, electrolysis or any other filler sculpting treatment within the last 2 weeks).
- Substance of the following statement: "This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety."
- Do not use this heater in the immediate surroundings of a bath, a shower or a swimming pool.
- Do not use this heater if it has been dropped.
- Do not use if there are visible signs of damage to the heater use this heater on a horizontal and stable surface.
- WEEE disposal information

Contraindications and other warnings

- If you have medical concerns, please consult a dermatologist first before use.
- Do not use the device on severe nodular or cystic acne.
- Do not use the device on moles, warts, ingrown hairs or irritated skin.
- Do not use the device if you are sensitive to light.
- Do not use the device if you have had sunburn or have burnt your skin within the last two weeks.
- Do not use the device directly on your eyes.
- Do not use the device more than 3 times in one day on the same area. Excessive use may burn or damage the skin. Ensure at least four hours between treatments.
- If you experience any skin discomfort after using the device (skin feels tight, uncomfortable, irritated, red, burning, etc.), notice unexpected side effects or prolonged redness, discontinue use immediately and consult a doctor. However, minor temporary redness is normal.
- Do not look directly into the light source of the probe for long periods of time as this may cause eye damage.
- This device is not a toy and should not be used by children.
- Do not share your device with others as this may lead to cross-contamination of bacteria.
- Do not use the instrument if it is visibly damaged and do not attempt to dismantle or repair the instrument.

Product structure

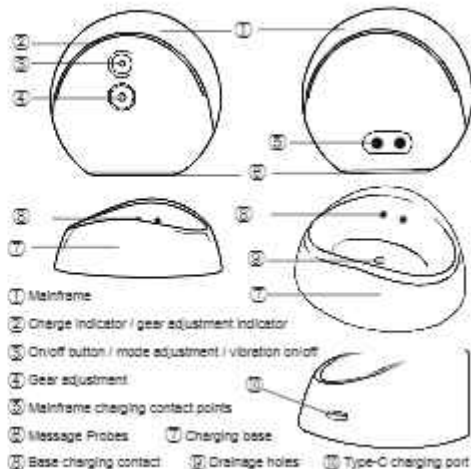

Functional Guidelines

1. The instrument has three modes: Stamp mode, Slide mode and EMS mode.
 2. After the instrument has been switched on by pressing and holding the on button, you can select the desired function mode by mode switching.
 3. After the instrument is switched on, you can select the intensity that suits you by the gear switching key.
 4. Press and hold the on/off key for 3 seconds to switch on or off the instrument.
- *Automatic protection shutdown function: When you do not press any buttons, the instrument will automatically switch off after 5 minutes of continuous use.

5. How do I turn on vibration?

No matter which mode is selected, the vibration will not be turned on, if you want to turn on the vibration, please double click the switch button; if you double click the switch button again, the vibration will be cancelled; switching modes will also turn off the vibration automatically.

Explanatory note on indicator lights

Charging instructions


1-8 indicator area lit

Gear indication

Low gear: 3-4 light

Mid-range: 2-5 bright

High level: 1-8 fully lit

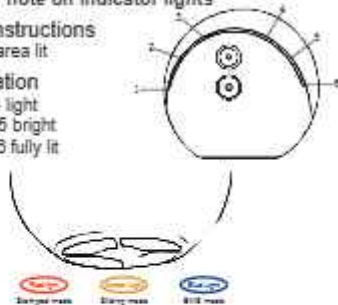

Functional Guidelines

Stamped mode

1. Principle and efficacy

This mode RF radio frequency will work simultaneously with the nano red light.

RF waves enter the subcutaneous tissue, causing the natural resistance movement of the subcutaneous tissue to generate heat energy, using the true collagen fibres in the dermis will

produce contraction, which has the dual effect of tightening the skin and promoting collagen regeneration. This has the dual effect of tightening the skin and promoting collagen regeneration, as the immediate contraction of collagen stimulates the dermis to secrete more new collagen fibres for regeneration. The immediate contraction of collagen stimulates the secretion of new collagen fibres in the dermis to regenerate, which in turn lifts the fibrous tissue of the skin and increases the thickness and density of the dermis, filling in wrinkles and improving the appearance of the skin. This increases the thickness and density of the dermis, filling in wrinkles and improving sagging.



LED nano red light uses a cosmetic light source that improves the microcirculation of the lymphatic system and stimulates intracellular mitochondrial. The specific wavelength of led red light can neutralize the free radicals produced in the cells and stimulate fibroblasts to produce. The specific wavelengths of led red light neutralise the free radicals produced in the cells, stimulate collagen production in the fibroblasts and accelerate blood circulation to achieve the effect of photorejuvenation.

The continuous vibrating waves generated by the instrument massage the skin, tighten it, promote nutrient absorption and effectively prevent wrinkles and sagging skin. The continuous vibrations produced by the device massage the skin, tighten it, promote nutrient absorption and effectively prevent wrinkles and sagging skin.

2 How to use

Step 1: Press and hold the on button to turn on the machine, the default start stamp mode and the default is 1 speed. According to your own skin the level of stimulation received, adjust the gear. The higher the setting, the better the effect.

Step 2: After applying the nutritional essence (or when applying the mask), use the massage probe to contact the skin and perform the stamped press Pressing operation.

Step 3: Press the seal in order on the face. When 8 seconds have elapsed, a beep will sound, indicating that the area has been stamped, change to the next stamped area until the whole face has been stamped, ending the stamped beauty round.

Description: The RF waves generated by this mode will automatically drop after the "beep" tone and enter a semi-hibernation. After the "beep", the RF energy will automatically drop and enter a semi-sleep state. Therefore, even if the instrument is not removed from the skin in time after the "beep" reminder, no damage will be caused to the skin.

Sliding mode

1.Principle and efficacy

RF radio frequency will work simultaneously with nano-amber cosmetic light.

RF waves produce heat energy that increases the thickness and density of the skin's dermis, repairing the reticular fibres, fading wrinkles and improve sagging.



LED Nano Amber Light uses 588-598nm wavelength to improve the microcirculation of the skin, increase skin elasticity, soothe wrinkles, stretch marks and slow down skin ageing. Warning: The RF waves generated by this mode will continuously heat the dermis of the skin, so please use the massage probe gliding operation all the way through when it is pressed against the skin, without staying too long, for a maximum of 8 seconds, do not stay in the Do not leave the probe on the same area of the skin for too long to avoid burning the skin.

2.How to use

Step 1: Press and hold the on button to turn on the device, the default is to start the stamp mode (red light) first, tap the on button, that is to cut switch to [slide mode] (amber light), the default is 1 gear, you can according to their own skin to accept the degree of stimulation, adjust the gears; the higher the gear the better the effect.

Step 2: After applying the nutritional essence (or when applying the mask), use the massage probe to contact the skin and perform a slow sliding operation. After applying the nourishing essence (or while applying the mask), use the massage probe to touch the skin and perform a slow gliding action.

EMS mode

1.Principle and efficacy

The EMS microcurrent will work simultaneously with the nano blue beauty light.

The EMS microcurrent revitalises the facial muscles and epidermal nerves, encouraging a more effective contraction of the facial skin. The skin's puffiness is reduced and the skin is lifted. The LED nano blue light uses a 465nm light source.

The LED nano blue light with 465nm wavelength regulates the oil-water balance of the skin and has a cleansing and sterilising effect. It has a cleansing and sterilising effect on the skin and eliminates acne.

2.How to use

Step 1: Press and hold the power on button to turn on the machine, the default start stamp mode (red light), tap the on/off button to switch to blue light, that is EMS mode, the default is 1 gear, according to their own skin to accept the degree of stimulation, the make the gear adjustment, the larger the gear, the better the beauty effect.

Step 2: After applying the nutritional essence (or when applying the mask), use the massage probe to touch the skin and perform Slow glide operation.

Care advice

Gear recommendations: can be tailored to individual circumstances.



*We recommend [Stamp Mode] for 5 minutes to promote collagen regeneration and [Slide Mode] with the mask. Use the remaining essence or gel for 5 minutes to lift and firm.



Stamped Anti-Aging
mode in 5 minutes



Sliding face sculpting
mode for 5 minutes

[Stamped Mode]:

Recommended for use 1-4 times a week depending on individual skin type.

- Normal skin type - 3 times a week for 10-15 minutes
- Dry skin - 2 times a week for 10-15 minutes
- Oily skin - 15-20 minutes 3-4 times a week
- Sensitive skin - 1 time per week for 5-10 minutes

[Slide Mode]:

•It is recommended to use the device once every 2 days. The [Stamped Anti-Aging Mode] and [Sliding Face Sculpting Mode] can be used together for a comprehensive beauty treatment. It is recommended to use the [Stamped Anti-Aging Mode] for 5 minutes to promote collagen regeneration, followed by the [Sliding Face Sculpting Mode] for 5 minutes using the remaining essence of the mask to lift and firm the skin. Alternatively, the [Sliding Mode] can be used with the remaining essence or gel of the mask for 5 minutes to achieve a lifting and firming effect.

[EMS mode]:

•It is recommended to use the device once every 2 days. The [EMS mode] enhances the beauty treatment and can be used in conjunction with the first two modes or on its own, in any order, to accelerate skin cell metabolism, inhibit bacteria, and help eliminate acne.

For optimal results, the recommended order of application for all three modes on the face is: 1. Left cheek—2. Right cheek 3. Forehead

Care steps

Step 1: Makeup removal and cleansing

- Remove all make-up and cleanse the face and neck thoroughly to ensure the skin is full clean.



Step 2: Apply a gel or mask

[Stamped Anti-Aging Mode]:



- Remove the mask, unfold it and apply it to the area of the face where it is to be used (the mask should fit flat and snugly on the face).
- or apply a gel or mask to the face or apply directly onto the probe.

Step 3: Facial

- Turn on the device and select the desired mode. You can use each of the 3 modes once per session.
- If no key operation is performed during use, the device will automatically switch off after 5 minutes of inactivity. To continue, simply turn the device back on and resume use.



Step 4: Cleaning and hydration

- At the end of the treatment, wash your face with water and apply the appropriate amount of skin care products to hydrate and moisturise according to your personal needs.


Step 5: Machine maintenance

- At the end of the treatment, wipe any residue from the beauty device with a slightly damp soft cloth, dry it and place back on the charging stand.
- If there is a lot of residue, you can rinse it off directly with water, rinse it off, dry it with a dry cloth and put it back on the charging stand. A small amount of water residue in the device crevices is normal.
- The charging cradle is not waterproof and cannot be rinsed with water.
- Can be rinsed
However, it must be dried after rinsing before use. Otherwise there will be residual moisture which will affect the experience.



How to charge

1. Connect the USB charging cable to a USB charging head with an output of 5V, (you can also access a computer or mobile power USB). Insert the Type-C port into the charging base connector.
2. Align the metal contacts on the back of the main unit with the metal contacts of the charging cradle and tilt it into the charging base. The arc indicator on the front of the main unit will light up in white and will be breathing, indicating the start of charging.
3. When charging, the arc-shaped charging indicator light will be breathing and undulating, when fully charged, the arc-shaped light will be on for a long time and will no longer be undulating and flashing. When fully charged, the arc light will remain on for a long time and will no longer fluctuate and flash.
4. When fully charged, if no operation is performed on the device, the indicator light will automatically turn off after approximately 30 minutes, and the device will switch-off to prevent interference. If the device is picked up again, the indicator will remain on for approximately 30 minutes before turning off, even if no operation has been performed.

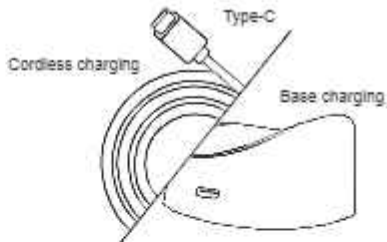

Caution:

Before charging, make sure that the charging slot is dry and free of water and skin care products. Do not charge the instrument in a place where it may be immersed in water. Do not charge the instrument in a place where it may be submerged in water. Do not turn on the switch while the instrument is charging.

Cleaning and maintenance of instruments

- Do not use wire balls, abrasive cleaners or corrosive liquids (e.g. acetone) to clean the instrument.
- Clean the instrument thoroughly with water or alcohol after each use and allow to air dry to ensure performance.

Comment question and solution:

Questions	Speculative reasons	Solutions
Unable to start the device	Battery depletion	Charging is complete before use
	Battery life exhaustion	Contact the seller for repair or replacement
Automatic equipment shutdown	No operation for more than 5 minutes	Rebooting
	Battery depletion	Charging is complete before use
	Insufficient amount of gel	Replenish the amount of gel on the face
Unable to charge	Using a mismatched Charging equipment	Use a charging plug with 5V output (or rechargeable battery)
Feeling too hot on the face	Over-gearred	Reduced gear strength
	Beauty device in local stays still	Slightly faster product movement
Feeling that your face is not hot enough	Gear too low	Upgraded gear strength
	Beauty device sliding too fast	Slightly reduced product movement speed
	Insufficient amount of gel	Replenish the amount of gel on the face
	Not enough energy storage when the unit is first switched on	Let the product run a little longer

Questions	Speculative reasons	Solutions
Current tingling sensation	EMS gears are too strong	Downshift
	Device does not fit snugly against the skin	Hold the device close to the skin
Sudden feverish tingling sensation	Only the edge of the probe is pressed against the skin	Place all 4 electrodes flat against the skin at the same time
Skin after application Slightly red	RF wave action of Normal phenomenon	It may return to normal in the short term, if it does not recover or have other adverse reactions, please seek medical attention
When in use water flows out	Product not drained of water	Pat the product on a dry cloth, drain or leave it to dry in a ventilated area. Ensure fully dry before using
In use there is "beep" beep	Different models	Only the stamping mode has a "beep" beep, the rest of the modes do not

Disclaimers

The use of this product is at the user's own risk. No liability or responsibility is accepted for any physical or non-physical damage, injury, or harm caused directly or indirectly by the use of the product. Additionally, the right is reserved to modify these instructions and change the contents at any time, without obligation to notify any individual of such modifications or changes.

Disposal methods

Please dispose of the product in accordance with local regulations for waste separation and recycling. Do not dispose of electronic equipment in regular household waste. Instead, take the product to an appropriate collection point for recycling or disposal of electronic waste. In the EU, this product is subject to the Waste Electrical and Electronic Equipment (WEEE) Directive, which ensures that discarded electronic products are properly recycled to minimize environmental impact.

For more information on recycling points, please contact your local waste management authority.



Product specifications

Product Model: 1607

Product power: $\leq 10W$

Frequency: 1000KHz(1MHz)

Battery capacity: 1100mA

Charging time: 3H

Input voltage: 5V

USB charging port: Type-C

Waterproof rating: IPX8

- Stamped mode battery working time about 60 minutes
- Sliding mode: battery working time about 60 minutes
- EMS mode: battery working time about 120 minutes



Designed & Distributed by:
JINJIANG CO., LTD.
Room: 1606, Zhenyuan Plaza, Zhongyuan
Nanjing, 210000, China
T: 025-26140000
Email: jiangjiang@jtc.com.cn

MADE IN CHINA

